

## Bowie and Dick Test

- Conforms to EN ISO 11140-4
- Independently validated to mirror the performance of the original Towel pack
- Diagnostic capabilities:
  - air removal
  - wet steam
  - super heated steam or presence of non-condensable gases
- Clear and accurate colour change
- Easy to interpret – easy to record
- Non-toxic, single use, disposable
- Lead free



### Bowie and Dick Daily Steriliser Test (B&D Test) The ParaSure B&D test pack fulfils the requirements of both EN 285 and EN ISO 17665

- The test is performed at the beginning of each working shift in every vacuum assisted steam steriliser.
- This non-toxic test is not a test of sterility but a test of the mechanical function of the steriliser to ensure rapid and even steam penetration into a small porous load.

Order Details

**PS1020**

Bowie and Dick Test pack  
Box of 20


## Daily Bowie & Dick Test Pack



Wrapper



Reverse of Indicator Sheet

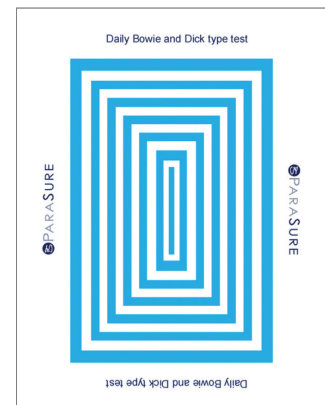

Unprocessed Test Pack Indicator Sheet

### TECHNICAL SPECIFICATIONS

<b>Name:</b>	Daily Bowie & Dick Test Pack, steam penetration
<b>Classification:</b>	EN-ISO 11140-1 (Class 2), EN-ISO 11140-4
<b>Operating temp:</b>	134°C - 137°C
<b>Time:</b>	3.3 min - 3.5 min
<b>Storage conditions:</b>	10°C - 30°C
<b>RH:</b>	30 - 70%
<b>Lot no.</b>	yes
<b>Prod date:</b>	yes
<b>Exp date:</b>	3 years
<b>Correct symbols:</b>	in accordance with EN-ISO 11140-4
<b>Indicator:</b>	both internal and external, external = process indicator, Internal = diagnostic indicator
<b>Disp/single use:</b>	yes
<b>Information on pack:</b>	
Cover	wrapper
Filtration medium	barrier system
<b>Extra information:</b>	
Sheet	silk screen printed
Coated	yes
Non-toxic/Lead free	yes
Lot no.	yes
<b>Recordable Information:</b>	
Hospital	yes
Department	yes
Autoclave no	yes
Cycle no	yes
Date	yes
Operator	yes
Supervisor	yes
Pass/fail result	yes
Colour change	blue to black

Address: Squirrels | Kithurst Park | Storrington | West Sussex | RH20 4JH | England |  
Tel: +44 (0) 1903 741390

Warehouse Address: Units B10/B11  
Oakendene Estate | Bolney Road | Cowfold | West Sussex | RH13 8AZ | England