




PARASURE

Bowie and Dick Test

- Conforms to EN ISO 11140-4
- CE Mark
- GTIN numbered
- Independently validated to mirror the performance of the original Towel pack
- Diagnostic capabilities:
 - air removal
 - wet steam
 - super heated steam or presence of non-condensable gases
- Clear and accurate colour change
- Easy to interpret – easy to record
- Non-toxic, single use, disposable
- Lead free

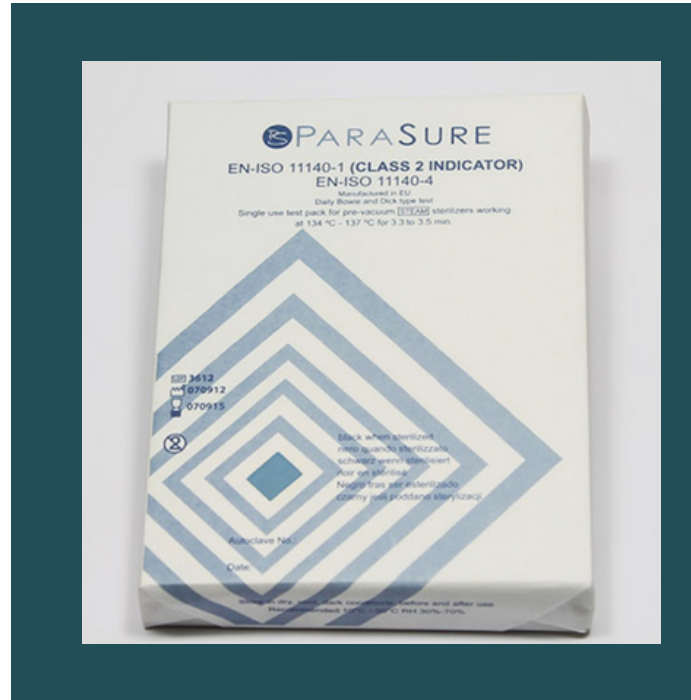

Bowie and Dick Daily Steriliser Test (B&D Test) The ParaSure B&D test pack fulfils the requirements of both EN 285 and EN ISO 17665

- The test is performed at the beginning of each working shift in every vacuum assisted steam steriliser.
- This non-toxic test is not a test of sterility but a test of the mechanical function of the steriliser to ensure rapid and even steam penetration into a small porous load.

Order Details

PS1020

Bowie and Dick Test pack
Box of 30

Registered Address:

ParaSure Limited Units B10/11 Oakdene
Estate Bolney Rd Cowfold West Sussex RH13
8AZ England

+44 (0)1403 864777

Ref: ParaSure Tech Sheets 02/13 Rev.D



WWW.PARASURE.CO.UK



WWW.PARASUREORGANICS.COM



PARASURE

Daily Bowie & Dick Test Pack

TECHNICAL SPECIFICATIONS

Name	Daily Bowie & Dick Test Pack, steam penetration
Classification:	EN-ISO 11140-1 (Type 2), EN-ISO 11140-4
Operating temp:	134°C - 137°C
Time:	3.3 min - 3.5 min
Storage conditions:	10°C - 30°C
RH:	30 - 70%
Lot no.	yes
Prod date:	yes
Exp date:	3 years
Correct symbols:	in accordance with EN-ISO 11140-4
Indicator:	both internal and external, external = process indicator, internal = diagnostic indicator
Disp/single use:	yes
Information on pack:	
Cover	Wrapper
Filtration medium	barrier system
Extra information:	
Sheet	silk screen printed
Coated	yes
Non-toxic/Lead free	yes
Lot no.	yes
Recordable Information:	
Hospital	yes
Department	yes
Autoclave no	yes
Cycle no	yes
Date	yes
Operator	yes
Supervisor	yes
Pass/fail result	yes
Colour change	pink to black

PARASURE

SITE KRANKENHAUS / OSPEDALE / HOSPITAL

AUTOClave NO / STERILISATOR NO / AUTOCLAVE NO / Nº ESTERILIZADOR

DEPARTMENT ABTEILUNG / REPARTO / DEPARTAMENTO

CYCLE NUMBER / CYCLOS NUMBER / NUMERO DI CICLO / Nº CICLO

SUPERVISOR LEITER / SUPERVISOR / SUPERVISOR

OPERATOR AUSGEFÜHRT VON / OPERATORE / OPERARIO

DATE DATUM / DATA / FECHA

RESULTS ERGEBNIS / RISULTATO / RESULTADO

Wrapper



Reverse of Indicator Sheet



Unprocessed Test Pack Indicator Sheet

For further information about any of our products, please visit:
www.parasure.co.uk



WWW.PARASURE.CO.UK

Registered Address:

ParaSure Limited Units B10/11 Oakendene
 Estate Bolney Rd Cowfold West Sussex RH13
 8AZ England

+44 (0)1403 864777



WWW.PARASUREORGANICS.COM